


Rhino Rocket[®]

Nasal Packing

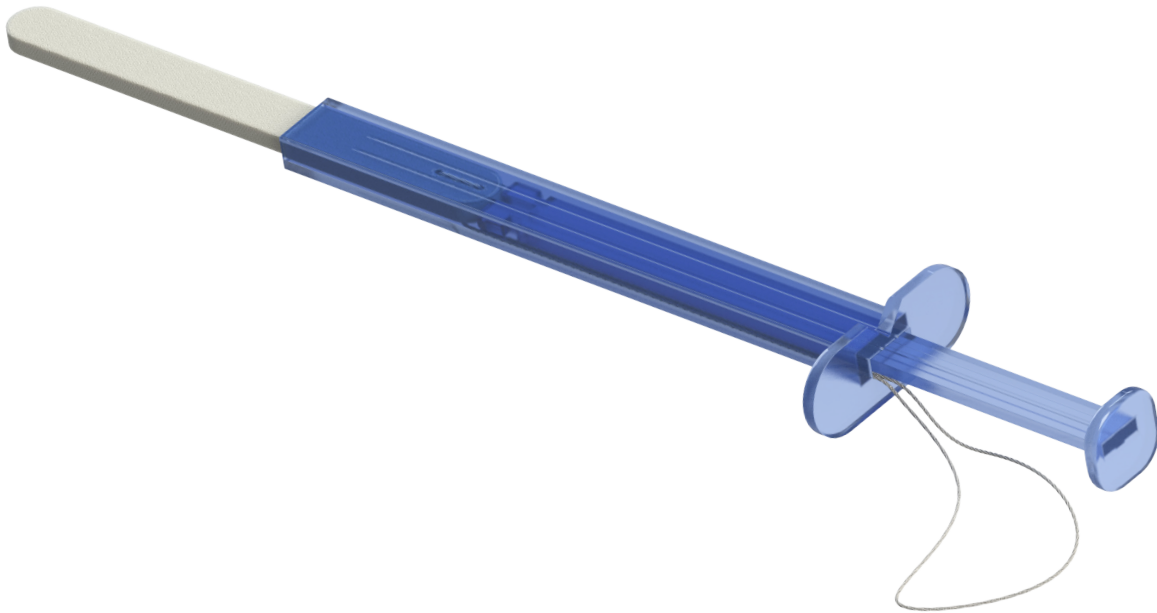

KEY PRODUCT FEATURES



- ✓ PVA Expandacell[®] foam for epistaxis management after trauma or surgery with its patented applicator
- ✓ Protects foam from moisture prior to insertions
- ✓ Aids in placement and insertion
- ✓ Minimizes blood containment concerns
- ✓ 30 seconds to insert
- ✓ Provides gentle, but firm mucosal compression on the tamponade of anterior or middle segment epistaxis
- ✓ When in contact with moisture, foam absorbs and swells to six times its compressed diameter
- ✓ String attached for retention and easy removal

- ✓ Eight sizes to treat pediatric to posterior bleeds
- ✓ Sterile

Rhino Rocket®

Nasal Packing

Rhino Rocket® Nasal Packing			
Slimline Rhino Rocket®			
Code	Description	Dimensions cm (LxWxH)*	QTY/Box
11S-S0300-08AS	Small	3 x 1 x 2	8
11S-S0500-08AS	Medium	5 x 1 x 2	8
11S-S0800-08AS	Large	8 x 1 x 2	8
11S-S1000-08AS	Extra Large	10 x 1 x 2	8
Wideline Rhino Rocket®			
Code	Description	Dimensions cm (LxWxH)*	QTY/Box
11W-W0300-08AS	Small	3 x 1.5 x 2	8
11W-W0500-08AS	Medium	5 x 1.5 x 2	8
11W-W0800-08AS	Large	8 x 1.5 x 2	8
11W-W1000-08AS	Extra Large	10 x 1.5 x 2	8
Rhino Rocket® Starter Pack			
Code	Description	Dimensions cm (LxWxH)*	QTY/Box
11-STARTERRHINO	One of each size	3 x 1 x 2 5 x 1 x 2 8 x 1 x 2 10 x 1 x 2 3 x 1.5 x 2 5 x 1.5 x 2 8 x 1.5 x 2 10 x 1.5 x 2	8

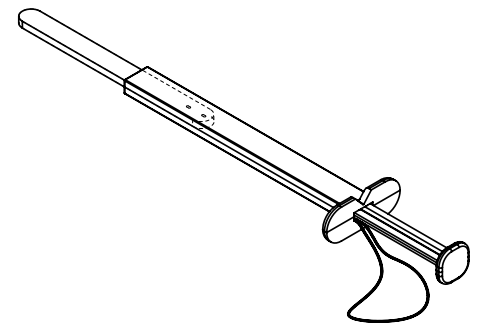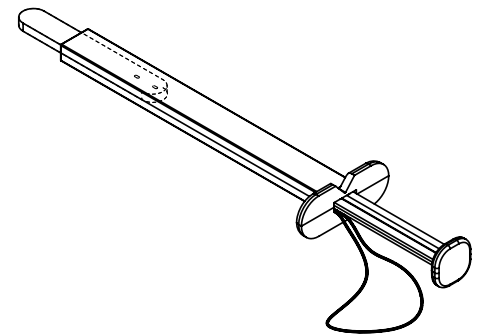

*Dimensions listed are for hydrated state and for insertion position