



# AbsorbENT™

## Nasal & Sinus Packs

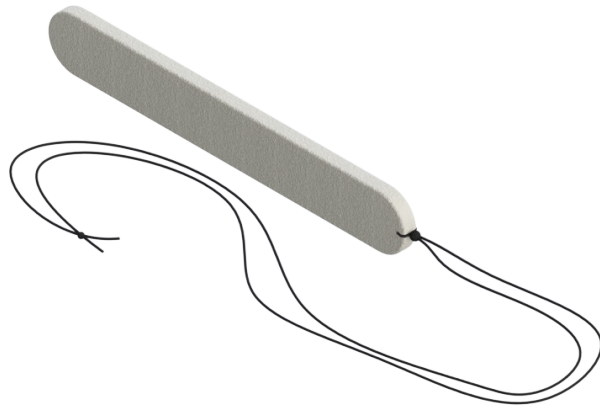

Compressed

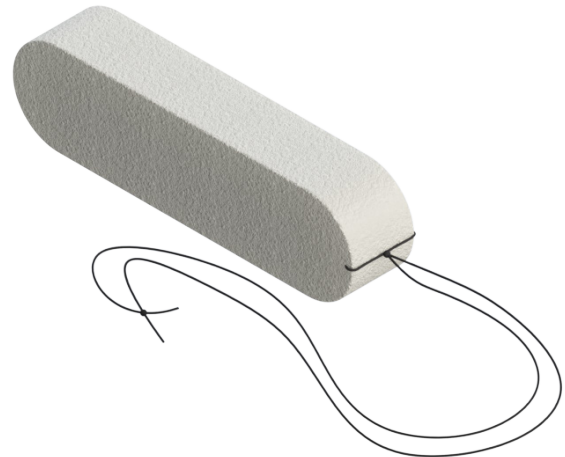


Hydrated

### KEY PRODUCT FEATURES

- ✓ Nasal and sinus packs utilize three-dimensional interconnected open cells that are soft and pliable for epistaxis management
- ✓ Placed into the operative area and when in contact with blood or other body fluids, they absorb the fluid and aid in the control of bleeding
- ✓ Provides gentle, but firm mucosal compression on the tamponade of anterior or middle segment epistaxis
- ✓ Can prevent lateralization of the turbinate
- ✓ Sterile Polyvinyl Alcohol (PVA) foam
- ✓ When in contact with moisture, foam absorbs and swells to six times its compressed diameter
- ✓ String attached for retention and easy removal
- ✓ Soft, lint and fiber free
- ✓ Rapid wicking properties

- ✓ Even and fast expansion
- ✓ Single-use

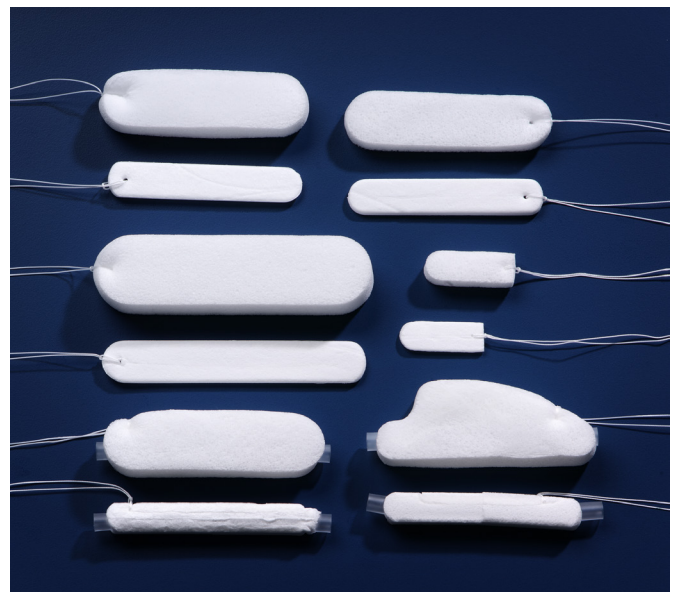

# AbsorbENT™

## Nasal & Sinus Packs


### AbsorbENT™ Nasal and Sinus Packing

#### Nasal Packs


Code	Description	Dimensions cm (LxWxH)*	QTY/Box
<b>RH-2315-20</b>	Small	5.5 x 1.5 x 2	20
<b>RH-7401-10</b>	Small	4.5 x 1.5 x 2	10
<b>RH-7408-10</b>	Small with Airway	4.5 x 1.5 x 2	10
<b>RH-7404-10</b>	Large	8 x 1.5 x 2	10
<b>RH-7405-10</b>	Large with Airway	8 x 1.5 x 2	10
<b>RH-7409-10</b>	Large Anatomical	8 x 1.5 x 3	10
<b>RH-7410-10</b>	Large Anatomical with Airway	8 x 1.5 x 3	10
<b>RH-7411-10</b>	Large Anatomical - Thin	8 x 1 x 3	10
<b>RH-7406-10</b>	Large - Thin	8 x 1 x 2	10
<b>RH-2310-10</b>	Extra Large	10 x 1.5 x 2.5	10

#### Sinus Packs

Code	Description	Dimensions cm (LxWxH)*	QTY/Box
<b>RH-7530-10</b>	Sinus .6W. Post-op	3.5 x 1.2 x 0.6	2/pack, 10 packs/box
<b>RH-7532-10</b>	Sinus .9W. Post-op	3.5 x 1.2 x 0.9	2/pack, 10 packs/box

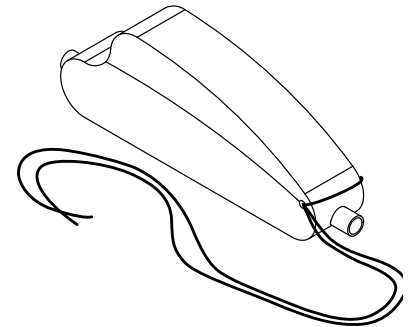

Anatomical with Airway

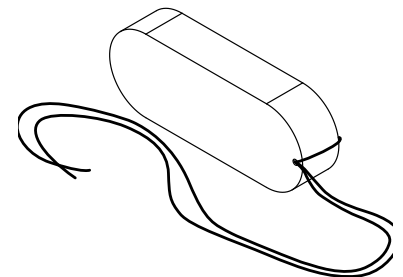

Sinus

\*Dimensions listed are for a hydrated state.

